



Canadian Friesian Horse Association

Box 2032 Athabasca, AB., T9S 2B6
Phone: 780 675 5927 Fax: 866 594 3120
Email: info@canadianfriesianhorse.ca
www.canadianfriesianhorse.ca

Application for Microchip Implant

This form must be completed when a horse is micro-chipped at any time other than at a Keuring. 'AVID and TROVAN' microchips should be used. If your vet can't supply the chip, please contact us and we will mail the microchip to you at a cost of \$20.00.

Breed of Horse - Friesian ____ Arabo-Friesian ____ Friesian Sport-horse ____ Other ____

Name of Owner: _____ Phone #: _____

Address: _____ Postal Code: _____

Name of Veterinarian : _____

Phone #: _____

Address: _____

Postal Code: _____

Date of Implantation: _____

Signature of Veterinarian:

Horse to be Implanted

Name: _____ Registry: _____

Reg. #: _____

Gender: _____

Date of Birth: _____

Implanted Microchip number:

Place Sticker here:

Name: (please print) _____

Date: _____
(Owner, Lessee or Authorized Agent)

Signature: _____
(Owner, Lessee or Authorized Agent)

Please provide a copy of the instructions to the veterinarian who is doing the implantation.

A Certified veterinarian MUST implant microchip.

After implantation, the chip should be scanned and verified. Please make sure to record the number and place the sticker on this Microchip application. Save a copy of this form for your records.

Please send this original completed application with the fee of \$20.00 to CFHA head office. Address above.

Instructions to implant the Microchip

1. The microchip must be implanted in the left side of neck of the horse in the triangle of:

- Left hand holding on the crest, fingertips pointing toward the croup of the horse.
- Right hand holding left from the shoulder, fingertips pointing toward the crest of the horse.
- In the triangle where the fingertips of both hands are coming together is where the microchip should be implanted.

2. Microchips and DNA Kits can be purchased from the CFHA head office.

3. For detailed instructions go to the following webpage:

www.avidequineid.com Click right on top on: [info Library](#).

There is a video to download and photos from a veterinarian, on how to implant the microchip.



ADMINISTRATION

The injection procedure requires three(3) scans.

1st Scan the patient prior to the procedure to ascertain if the animal is already identified.

2nd Scan the AVID Chip prior to the procedure to verify the number.

3rd Scan the patient after the procedure to again verify and accurately record the number.

CANINE/ELINE PROCEDURE

1. The injection site is dorsal midline just caudal to a line drawn between the scapula. The chip should be deposited in the deep subcutaneous plane as close to, or in, the trapezius muscle, over the dorsal spinous process of the fifth or sixth thoracic vertebra.
2. Prepare the site for an aseptic injection.
3. "Tent" the skin. Place the needle at the base of the "tent" at a 45 degree angle from the dorsal midline.
4. Insert the needle fully until the hub comes in contact with the skin. Slowly depress the syringe plunger until it stops.
5. Apply pressure on the injection site as the needle is withdrawn.

AVIAN PROCEDURE

Your assistant should restrain the bird with one hand holding the head and neck and the other hand holding the feet and wing tips. It is best to have your assistant transfer the bird to two-hand restraint. The birds body should be held close to your handler to further immobilize the wings. IF YOUR VETERINARY ASSISTANT IS NOT EXPERIENCED IN PROPER AVIAN RESTRAINT, DO NOT ATTEMPT THIS PROCEDURE.


1. The injection site should be in the thickest portion of the breast muscle. The feathers overlying this area should be parted and wiped with alcohol.
2. Place the needle, bevel up at approximately a 50 degree angle from the birds body.
3. Insert the needle ONLY UNTIL THE BEVEL OF THE NEEDLE IS COVERED BY TISSUE (or mostly covered in smaller birds). Slowly depress the syringe plunger until it stops.
4. Apply pressure on the injection site as the needle is withdrawn.

EQUINE PROCEDURE

1. The injection site is the left side of the neck, half the distance between the poll and withers, one inch below the midline in the nuchal ligament.
2. Shave the injection area and prepare the site for an aseptic injection.
3. Inject the microchip directly into the prepared site perpendicular to the neck. Slowly depress the syringe plunger until it stops.
4. Apply pressure on the injection site as the needle is withdrawn.

PLEASE READ AND FOLLOW INSTRUCTIONS

Fig 1

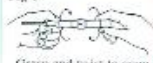

Grasp and twist to open.
CAUTION: DO NOT DEPRESS PLUNGER UNTIL NEEDLE IS INJECTED.

Fig 2



Grasp needle cover and slide outer sleeve down.

Fig 3



Grasp outer sleeve and remove the syringe assembly. DO NOT pull on plunger.

Fig 4



Insert the needle cover into the small side of the outer sleeve. Press firmly.

Fig 5



Remove SUIDS assembly from the needle cover.
CAUTION: Needle sharp!

DO NOT RECAP EXCEPT AS PERMITTED BY OSHA, CDC OR OTHER REGULATIONS. DISPOSE IN APPROVED MANNER ONLY.

AVID Identification Systems, Inc

3185 Hammer Avenue Norcross, California 92863 USA Tel: 41-951-371-7505 Fax: 41-951-737-8967 email: info@AvidID.com www.AvidID.com

USA 1-800-336-2843